



NOCCARC V310

V E N T I L A T O R

With
Integrated
HFNC

Easy to use

Reliable

Versatile



Noccarc V310 offers an advanced, compact yet reliable ICU ventilation solution. Carefully crafted with the combination of cutting edge research and advanced manufacturing technologies, the turbine based design enables highly efficient independent operation and a wide range of ventilation modes, including High Flow Nasal Cannula.

The IOT enabled system allows central data monitoring and logging remotely making **Noccarc V310** an accessible easy to manage option for the health professionals in the current case of Covid-19 pandemic.



Modes of Ventilation

Volume Controlled Ventilation	VC-CMV VC-AC VC-SIMV PRVC-CMV PRVC-AC PRVC-SIMV
Pressure Controlled Ventilation	PC-CMV PC-AC PC-SIMV
Spontaneous Breathing Support	PS/CPAP (PSV) with Apnea Backup Ventilation APRV HFNC
Therapy type	Invasive Ventilation Non-invasive Ventilation (NIV) O2 therapy
Types of Patients Supported	Adults and Paediatric
Apnea Backup Ventilation	
High Flow Nasal Oxygen Therapy	

Technical Specifications

Respiratory Rate (RR)	1- 60 BPM
Peak Inspiratory Pressure	60 cm H ₂ O
PEEP	0- 40 cm H ₂ O
Pressure Support	0- 50 cm H ₂ O
Inspiratory Time (Ti)	0.2- 5 s
PEEP	0- 50 cm H ₂ O
I:E Ratio	1:6 to 6:1
Tidal Volume (VT)	50- 2000 mL
Inspiratory Flow Rate	280 L/min (Internally restricted to 160 L/min)
Trigger Flow Sensitivity	1- 20 L/min
Volume Accuracy	2- 3% of full scale between 10 L/min- 80 L/m
O ₂ concentration (FiO ₂)	21- 100 %
FiO ₂ control	Closed loop electronic manipulation
Recruitment Tools	Inspiration Hold Expiration Hold
Leak Detection and Compensation	YES
SIGH	Programmable
O ₂ Therapy	Continuous flow 0- 100 L/min Concentration 21- 100 %
Nebulizer	for 5, 6, 7, 8, 9, 10 minutes Synced with Inspiratory Phase
Oxygen Input	45- 60 PSIG
Screen	10" TFT LCD Touchscreen
Dimensions Patient Unit	395 X 340 X 255 mm
Weight	~ 10.1 Kg Patient Unit ~ 25 Kg with trolley
Battery	Lithium, Internal
Battery Back-up Time	Upto 8 hours
Mains Input	200- 240V AC, 50- 60 Hz
Input Current	<2.00 A
Power Consumption	<150W

Compliant with IEC 60601-1 Standard

Monitored Parameters

Pressure Measurements	Peak Inspiratory Pressure (PIP) Positive end Expiratory Pressure (PEEP) Plateau Pressure (Pplat) Driving Pressure Mean Pressure (Pmean) Range: 0 to 100 cm H2O
Flow Measurements	Expiratory minute Volume (MVe) Inspiratory minute Volume (MVi) Spontaneous Expired minute Volume (Spont. MVe) Range 0 to 160 L/min
Volume Measurements	Inspiratory Tidal Volume (VTi) Expiratory Tidal Volume (VTe) Total Leak Volume
Other Measurements	Oxygen Concentration (FiO2) Respiratory Rate (RR) I:E Ratio Rapid Shallow Breathing Index (RSBI) Static Lung Compliance (Cstat) Dynamic Lung Compliance (Cdyn) Inspiratory Time (Ti) Expiratory Time (Te)
Trends	Storage of trends for last 48 hours of all parameters

Graphic Mode with Display

Flow vs Time

Pressure vs Time

Volume vs Time

Pressure-Flow Loop

Pressure-Volume Loop

Volume-Flow Loop





Alarms for Device and Patient Safety

Alarms for Patient Safety

PEEP (High / Low)
Respiratory Rate (High / Low)
Tidal Volume (High / Low)
FiO₂ (High / Low)
High Inspiratory Pressure
Patient Disconnection
Battery Low

Alarm Features

Color coded alarms based on priority
Adjustable Apnea Alarm time (Tapnea)
History of last 200 alarms

Unique Features




High Flow Oxygen Therapy HFNC Mode

Noccarc V310 HFNC mode when used with humidifier, is capable of delivering up to 100% humidified and heated oxygen at a flow rate of up to 100 LPM. High flow therapy increases alveoli recruitment and delivers high FiO_2 , while offering increased patient comfort.

UV-C Decontamination Chamber

The UV-C decontamination chamber mitigates the risk of transmission of airborne diseases by reducing the bacterial and viral load going into the ICU atmosphere from patient exhaled air.

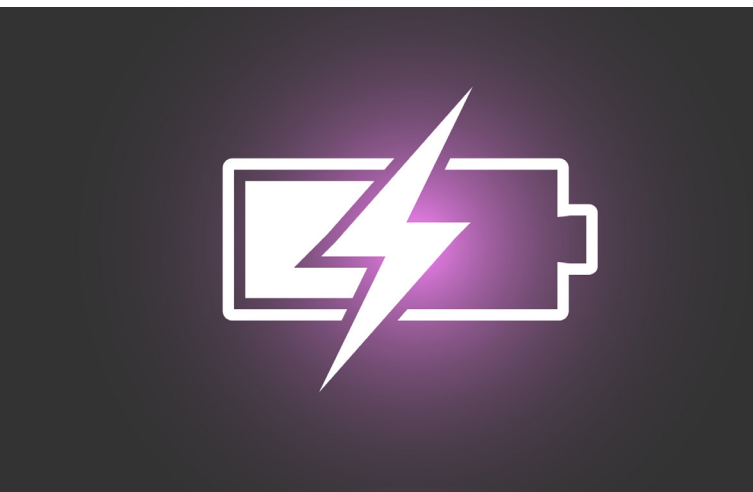

8 Hours Battery Backup

Noccarc V310 is a powerful turbine based ventilator with a battery backup of upto 8 hours. The compact, independent and power efficient design enables it to operate without compressed medical air or compressor, making it ideal for ICUs as well as transport.



Inbuilt Programmable Nebulizer Synced with Inspiratory Phase

Inbuilt programmable nebulizer provides quick and convenient way to administer medication to ventilated patients. The inspiratory phase synchronization algorithm reduces medication loss during expiration, and precisely delivers set tidal volume.

10 inch Touchscreen with User-friendly Interface

The integrated 10 inch Touchscreen Interface provides great monitoring and controlling capabilities. The customized interface is designed considering inputs from healthcare professionals and has been created for quick, easy and intuitive usage. A screen lock feature prevents accidental changes.





IoT Enabled Central Data Monitoring and Logging*

Central monitoring system incorporates innovative features such as flexible trending and display capabilities, advanced alarm management and cloud data logging. Combined with a user friendly web interface with seamless networking system via wireless LAN, Noccarc V310 offers safe and manageable control for healthcare professionals.

* Central monitoring system is not included with the ventilator and has to be ordered separately.

Testimonials

40, Sassoon Road, Pune - 411 001 (INDIA). Tel. : 020 - 6645 5100 (60 Lines), 2616 3391 (8 Lines)
Fax : 020 - 2616 4529 • E-mail : info@rubyhall.com • Website : www.rubyhall.com

DR. PRACHEE SATHE
M.D. FRCP, FCCM
Reg. No. 53396
Director & Head of the Dept. Critical Care Medicine
rubyicu2003@gmail.com
Phone No. : 020 66455243
www.rubyicu2003.com

Date: 26.06.2020


CERTIFICATE

This is to certify that Noccare V310, an ICU ventilator, developed by Nocca Robotics Pvt. Ltd. in collaboration with IIT Kanpur has been validated for clinical performance and functionality at Ruby Hall Clinic, Pune under the able guidance and supervision of Dr Prachee Sathe, Director-Dept. of Critical care Medicine and the team of intensivists in the ICU.



The device was subjected to the following tests:

1. Endurance Testing for 24 days
2. All the modes of operation, features and device settings as stated in the product brochure
3. Testing on a healthy volunteer
4. Testing on a critically ill patient suffering from pneumonia with respiratory failure

It is to further state that the ventilator's performance has been found satisfactory and Noccare V310 is good to be used as an ICU Ventilator.


Dr Prachee Sathe
 MD FRCP FCCM
 DIRECTOR ICU

PS-38-048

Medanta Institute of Critical Care and Anesthesiology


Dr. Yatin Mehta
MD, M.N.A.M.S, F.R.C.A, F.A.M.S, F.I.C.C.M, F.T.E.E
Chairman
Institute of Critical Care and Anesthesiology
yatin.mehta@medanta.org

Dated: 29th June 2020



TO WHOM IT MAY CONCERN

VENTILATOR NOCCARC V 310

The above ventilator was demonstrated to me on 20th June 2020. The NOCCARC V310 ventilator meets all the technical specifications and features as stated in the brochure. The ventilator performed satisfactorily and is fit to be used in the ICU.


Dr. Yatin Mehta
 M.D., M.N.A.M.S, F.R.C.A., F.A.M.S.
 Chairman
 Medanta Institute of Critical Care & Anesthesiology
 Medanta-The Medicity
 Sector-38, Gurgaon-122001
 Regn. No. DMAC-6690, HN-204499

Accredited by

Medanta - The Medicity
♦ Sector - 38, Gurgaon, Haryana, 122 001, India

Medanta - Mediclinic
♦ E - 18, Darkhosi Culture, New Delhi - 110 028
Tel. : 011-41419121/122 Fax : 011-41249314/33

Medanta - Mediclinic Cybercity
♦ 100 Feet, Building 90C, DE, Cyber City Phase II, Gurgaon 122 002
Tel. : 011-124-61441-422

Nocca Robotics Private Limited

Pune

Plot No-BGP 10/12, BSNL Chowk, Opposite Enzene
Biosciences, MIDC Bhosari, Pune- 411026

Delhi NCR

710, Logix City Centre, 7th Floor, Plot No. BW-
58, Sector- 32, Noida, Uttar Pradesh, 201301

Revision Date : 30/07/20